

Allergies and Hypersensitivities

Article Number: 52 | Rating: Unrated | Last Updated: Tue, Jul 19, 2016 at 8:22 AM

Overview

The form allows entering new allergy/hypersensitivity information and creating reports and reviews of all collected data.

Applicable Staff

The ID Nurse can enter allergies and hypersensitivities. All users have viewing access.

Usage

The form is client specific (non-episodic) and allows adding the information to active and inactive episodes or both on the first screen. Open the form from a hyperlink in the sidebar of one of the other forms that is connected to it, from within a client record, or the menu path Avatar CWS/Other Chart Entry/Allergies and Hypersensitivities. The form opens as shown below:

Avatar 2016

Home John D Preferences Lock Sign Out Switch Help SYSADM

JOHN DOE (000000013)
M, 60, 05/05/1955
Ht: 0' 1.0", Wt: 1 lbs, BMI: 703

Allergies (0)

Display shows a report of allergies. Choose whether to show results from the active episode, inactive, or both.

Allergies always show in the Client Data Bar. Review status shows when entered.

Update opens the data entry table to review, delete, or enter new info. Review info shows in the Client Data Bar above.

Include Allergies on Report: ☐ Active ☐ Inactive ☒ Both

Display

Allergies/Hypersensitivities Reviewed: ☐ Yes ☐ No

Update

Known Medication Allergies: ☐ Yes ☐ No

Submit

Online Documentation

AVPM (UAT) 04/27/2016 10:07 AM 100%

When the **Update** button is clicked, the following multiple iteration table opens to enter new data or review/delete previous entries. For detailed information regarding using multiple iteration tables, see this article. Using a Form

Avatar 2016

Home John D Preferences Lock Sign Out Switch Help SYSADM

JOHN DOE (000000013)
M, 60, 05/05/1955
Ht: 0' 1.0", Wt: 1 lbs, BMI: 703

Allergies (0)

Required field. Use T for today.

Date Recorded automatically populates with today's date

	Allergen/Reactant	Date Recorded	Status	Reactions	Reaction Severity	Onset	Date Recorded	Treatment	Comments
1	PEANUTS (MDX-2889) (SNOMED:256349002) (MDX-2889)	04/27/2016	Active (1)	Altered mental status&...	Moderate to sev...	Childhood (1)	04/27/2016	Epinephrine injection - EpiPen	Provide nut-free meals
2							04/27/2016		

Begin typing and press Enter. A pop up appears with matching results. Select and press TAB.

Status shortcuts:
1 - Active
2 - Inactive
3 - Resolved

Double click the box to show available results. Choose as many as appropriate and press TAB.

Reaction Severity search results:

Code	Description
SCT-399166001	Fatal
SCT-442452003	Life threatening severity
SCT-255604002	Mild
SCT-371923003	Mild to moderate
SCT-6736007	Moderate
SCT-371924009	Moderate to severe
NTST-NoReaction	No Reaction
SCT-24484000	Severe
NTST-Unknown	Unknown

Onset shortcuts:
1 - Childhood
2 - Adolescent
3 - Adult
4 - Late

Treatment


Treatment and Comments open a text box like this to enter information. Click OK when done to enter the info.

Click Save when complete to enter the info

New Row Delete Row Save Close/Cancel Copy/Paste Row Copy Cell Paste Cell

AVPM (UAT) 04/27/2016 10:38 AM

After clicking **Save** in the table above, the view returns to the Allergies and Hypersensitivities screen of the form. Click the blue **Submit** button to save the info into the record. When allergy/hypersensitivity data is entered, it will be shown in the **Client Data Bar** on all screens.

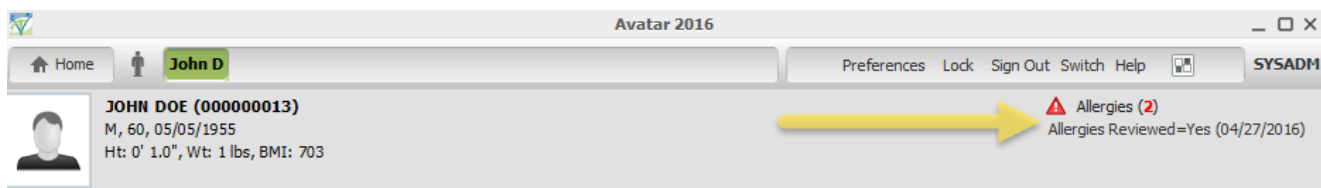

Tips & FAQ

Any useful tips, and questions provided by users after roll out.

Technical Information

Menu Location (Avatar CWS/Other Chart Entry/Allergies and Hypersensitivities)

SQL Table(s)

RAD Envelope (if applicable)

Related Forms and Reports (Individual Information Sheet, Medical History Summary)

Posted by: **John Cookson** - Tue, Apr 26, 2016 at 6:47 PM. This article has been viewed 3056 times.

Online URL: <https://www.mercykb.com/article.php?id=52>